

## Cytokines

### Human Recombinant IL-2 (CHO-expressed)

#### Interleukin 2

Catalog #	78036.1	10 µg
	78036	50 µg
	78036.3	100 µg
	78036.2	1000 µg



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

Interleukin 2 (IL-2) is a monomeric cytokine that was originally identified as a T cell growth factor (Gaffen & Liu). It binds to heterotrimeric receptors consisting of CD25, CD122, and CD132. Upon binding, it activates JAK3-, STAT5-, and AKT-dependent signaling pathways, which results in cellular proliferation and survival (Ma et al.). The majority of IL-2 is secreted by activated CD4+ and CD8+ T cells, although B cells and dendritic cells were found to produce IL-2 in small amounts. IL-2 downregulates immune responses to prevent autoimmunity during thymic development, influences the development of CD4+CD25+ regulatory T cells, and affects development of follicular helper T cells. IL-2 also controls inflammation by inhibiting Th17 differentiation (Banchereau et al.). High IL-2 levels in serum are associated with progression of scleroderma, rheumatoid arthritis, and gastric and non-small cell lung cancer, though no known disease can be directly attributed to the lack or excess of IL-2 (Gaffen & Liu).

## Product Information

Alternative Names:	Aldesleukin, Interleukin-2, T cell growth factor, TCGF
Accession Number:	P60568
Amino Acid Sequence:	APTSSSTKKT QLQLEHLLLD LQMILNGINN YKNPKLTRML TFKFYMPKKA TELKHLQCLE EELKPLEEVL NLAQSKNFHL RPRDLISNIN VIVLELKGSE TTFMCEYADE TATIVEFLNR WITFCQSIIS TLT
Predicted Molecular Mass:	15.3 kDa
Species:	Human
Formulation:	Lyophilized after dialysis against phosphate-buffered saline.
Source:	CHO


## Specifications

Activity:	The EC50 is $\leq 2$ ng/mL as determined by a cell proliferation assay using CTLL-2 cells. The specific activity is approximately $4.1 \times 10^4$ IU/µg as calibrated against the human recombinant IL-2 WHO International Standard (NIBSC code: 86/500).
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is $\leq 0.2$ EU/µg protein.

## Preparation and Storage

Storage:	Store at -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. As a general guide, do not store at 2 - 8°C for more than 1 week or at -20°C to -80°C for more than 2 months. Avoid repeated freeze-thaw cycles.

## Data

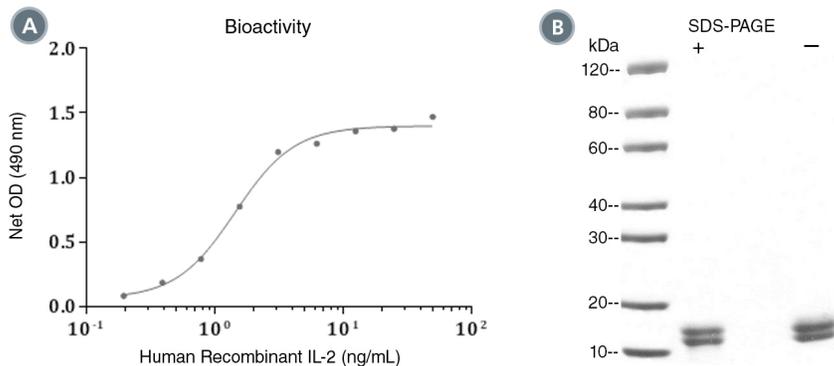

(A) The biological activity of Human Recombinant IL-2 (CHO-expressed) was tested by its ability to promote the proliferation of CTLL-2 cells. The EC50 is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC50 in the above example is 1.2 ng/mL.

(B) 2  $\mu$ g of Human Recombinant IL-2 (CHO-expressed) was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant IL-2 (CHO-expressed) has a predicted molecular mass of 15.3 kDa.

## Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit [www.stemcell.com/cytokines](http://www.stemcell.com/cytokines) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

- Banchereau J et al. (2012) From IL-2 to IL-37: the expanding spectrum of anti-inflammatory cytokines. *Nat Immunol* 13(10): 925–31.  
 Gaffen SL & Liu KD. (2004) Overview of interleukin-2 function, production and clinical applications. *Cytokine* 28(3): 109–23.  
 Ma A et al. (2006) Diverse functions of IL-2, IL-15, and IL-7 in lymphoid homeostasis. *Annu Rev Immunol* 24: 657–79.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2021 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.