

# EasySep™ Mouse Naïve CD8+ T Cell Isolation Kit

For processing 1 x 10<sup>9</sup> cells

Catalog #19858  
#19858RF RoboSep™

Negative Selection

Document #10000003738 | Version 02



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713  
INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM  
FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Description

Isolate untouched and highly purified naïve CD8+ T cells (CD3+CD8+CD44-CD62L+) from mouse splenocytes or lymph nodes by immunomagnetic negative selection. When using single-cell suspensions from other tissue types, this kit may require optimization.

- Fast, easy-to-use, and column-free
- Up to 98% purity
- Untouched, viable cells

This kit targets non-naïve CD8+ T cells for removal with biotinylated antibodies recognizing specific cell surface markers. Unwanted cells are labeled with biotinylated antibodies and streptavidin-coated magnetic particles and separated without columns using an EasySep™ magnet. Desired cells are simply poured off into a new tube. Isolated cells are immediately available for downstream applications, such as flow cytometry, culture, and cell-based experiments.

## Component Descriptions

COMPONENT NAME	COMPONENT #	QUANTITY	STORAGE	SHELF LIFE	FORMAT
EasySep™ Mouse Naïve CD8+ T Cell Isolation Cocktail	19858C	1 x 0.5 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS and 0.1% BSA.
EasySep™ Streptavidin RapidSpheres™ 50001	50001	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A suspension of magnetic particles in PBS.
EasySep™ Mouse FcR Blocker	18720	1 x 0.1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS, 0.1% BSA, and < 0.1% sodium azide.

BSA - bovine serum albumin; PBS - phosphate-buffered saline

Components may be shipped at room temperature (15 - 25°C) but should be stored as indicated above.

## Sample Preparation

For automated and standardized tissue processing, see STEMprep™ Tissue Dissociator (Catalog #100-2112) at [www.stemcell.com/stemprep](http://www.stemcell.com/stemprep). For manual processing, follow the steps below.

SPLEEN or LYMPH NODE

Disrupt spleen or lymph node in PBS or Hanks' Balanced Salt Solution (HBSS) containing 2% fetal bovine serum (FBS). Remove aggregates and debris by passing cell suspension through a 70 µm mesh nylon strainer (e.g. Catalog #27216). Centrifuge at 300 x g for 10 minutes and resuspend at 1 x 10<sup>8</sup> nucleated cells/mL in recommended medium.

Ammonium chloride treatment is not recommended when preparing the cells for separation.

## Recommended Medium

EasySep™ Buffer (Catalog #20144), RoboSep™ Buffer (Catalog #20104), or PBS containing 2% FBS and 1 mM EDTA. HBSS, Modified (Without Ca<sup>++</sup> and Mg<sup>++</sup>; Catalog #37250) can be used in place of PBS. Medium should be free of Ca<sup>++</sup>, Mg<sup>++</sup>, and biotin.

## Directions for Use – Manual EasySep™ Protocols

See page 1 for Sample Preparation and Recommended Medium. Refer to Tables 1 and 2 for detailed instructions regarding the EasySep™ procedure for each magnet.

**Table 1. EasySep™ Mouse Naïve CD8+ T Cell Isolation Kit Protocol**

		EASYSEP™ MAGNETS	
STEP	INSTRUCTIONS	 EasySep™ (Catalog #18000)	 “The Big Easy” (Catalog #18001)
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10 <sup>8</sup> cells/mL 0.25 - 2 mL	1 x 10 <sup>8</sup> cells/mL 0.5 - 8.5 mL
2	Add FcR blocker to sample.	10 µL/mL of sample	10 µL/mL of sample
3	Add sample to required tube.	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Catalog #38007)	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)
4	Add Isolation Cocktail to sample. NOTE: Do not vortex cocktail.	50 µL/mL of sample	50 µL/mL of sample
	Mix and incubate.	RT for 10 minutes	RT for 10 minutes
5	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds	30 seconds
6	Add RapidSpheres™ to sample.	100 µL/mL of sample	100 µL/mL of sample
	Mix and incubate.	RT for 5 minutes	RT for 5 minutes
7	Add recommended medium to top up the sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 2.5 mL	<ul style="list-style-type: none"> <li>• Top up to 5 mL for samples &lt; 4 mL</li> <li>• Top up to 10 mL for samples ≥ 4 mL</li> </ul>
	Place the tube (without lid) into the magnet and incubate.	RT for 2.5 minutes	RT for 2.5 minutes
8	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring the enriched cell suspension into a new tube.	Isolated cells are ready for use	Isolated cells are ready for use

RT - room temperature (15 - 25°C)

\* Leave the magnet and tube inverted for 2 - 3 seconds, then return upright. Do not shake or blot off any drops that may remain hanging from the mouth of the tube.

Table 2. EasySep™ Mouse Naïve CD8+ T Cell Isolation Kit Protocol

		EASYSEP™ MAGNET
STEP	INSTRUCTIONS	EasyPlate™ (Catalog #18102)
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10 <sup>8</sup> cells/mL 0.05 - 0.2 mL
2	Add FcR blocker to sample.	10 µL/mL of sample
3	Add sample to required tube.	Round-bottom, non-tissue culture-treated 96-well plate (e.g. Catalog #38018)
4	Add Isolation Cocktail to sample. NOTE: Do not vortex cocktail.	50 µL/mL of sample
	Mix and incubate.	RT for 10 minutes
5	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds
6	Add RapidSpheres™ to sample.	100 µL/mL of sample
	Mix and incubate.	RT for 5 minutes
7	Add recommended medium to top up sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 0.25 mL
	Place the plate (without lid) into the magnet and incubate.	RT for 2.5 minutes
8	Carefully pipette (do not pour) the enriched cell suspension into a new plate.	Isolated cells are ready for use

RT - room temperature (15 - 25°C)

## Directions for Use – Fully Automated RoboSep™ Protocol

See page 1 for Sample Preparation and Recommended Medium. Refer to Table 3 for detailed instructions regarding the RoboSep™ procedure.

**Table 3. RoboSep™ Mouse Naïve CD8+ T Cell Isolation Kit Protocol**

STEP	INSTRUCTIONS	RoboSep™ (Catalog #21000)	
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10 <sup>8</sup> cells/mL 1 - 8.5 mL NOTE: If starting with 0.5 - 1 x 10 <sup>8</sup> cells, resuspend cells in 1 mL.	
2	Add FcR blocker to sample.	10 µL/mL of sample	
3	Add sample to required tube.	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)	
4	Select protocol.	Mouse Naïve CD8+ T Cell Isolation 19858	
5	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds	
6	Load the carousel.	Follow on-screen prompts	
	Start the protocol.	Press the green "Run" button	
7	Unload the carousel when the run is complete.	Isolated cells are ready for use	

## Notes and Tips


### ASSESSING PURITY

For purity assessment of naïve CD8+ T cells (CD3+CD8+CD44-CD62L+) by flow cytometry, use the following fluorochrome-conjugated antibody clones:

- Anti-Mouse CD3e Antibody, Clone 145-2C11 (Catalog #60015), and
- Anti-Mouse CD8a Antibody, Clone 53-6.7 (Catalog #60023), and
- Anti-Mouse CD62L (L-Selectin) Antibody, Clone MEL-14 (Catalog #60109), and
- Anti-mouse CD44 (Ly-24) antibody, clone 5035-41.1D

The anti-mouse CD44 (Ly-24) clone 5035-41.1D is not blocked by the anti-CD44 clone used in EasySep™ Mouse Naïve CD8+ T Cell Isolation Cocktail. The 5035-41.1D clone only recognizes the Ly-24.2 isoform, which is expressed by C57BL/6, C57BL/10, C57/L, C58A, AKR, 129, SJL, NZB, C3H, CE, and CBA/H mouse strains. The anti-mouse/human CD44 clone IM7 is not recommended for assessing purity, as it is blocked by the anti-mouse CD44 antibody clone used in EasySep™ Mouse Naïve CD8+ T Cell Isolation Cocktail.

## Data

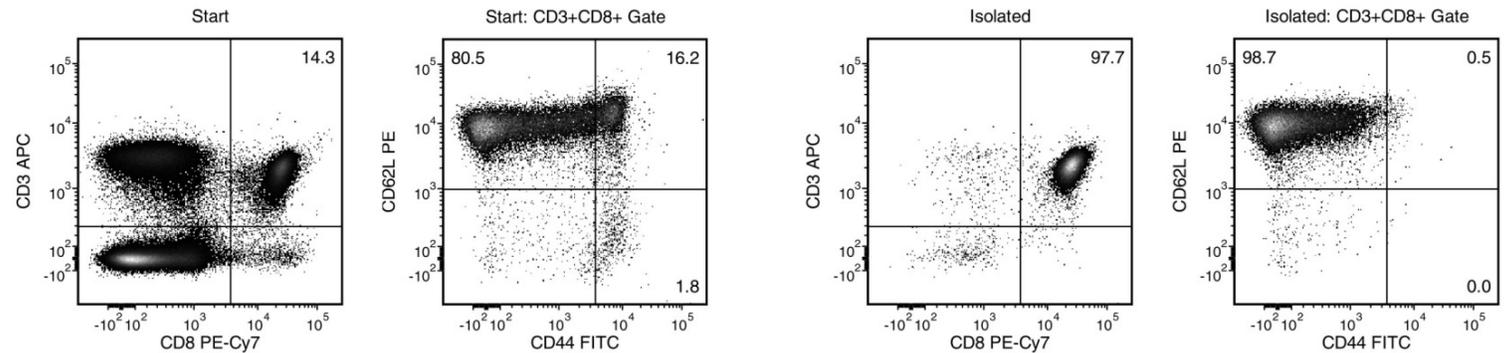

**Figure 1. Typical EasySep™ Mouse Naïve CD8+ T Cell Isolation Profile**

Starting with splenocytes from an uninfected mouse, the naïve CD8+ T cell content (CD3+CD8+CD44-CD62L+) of the isolated fraction typically ranges from 92 - 98%. In the above example, the purities of the start and final isolated fractions are 11.5% and 96.4%, respectively.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED. FOR ADDITIONAL INFORMATION ON QUALITY AT STEMCELL, REFER TO [WWW.STEMCELL.COM/COMPLIANCE](http://WWW.STEMCELL.COM/COMPLIANCE).

Copyright © 2025 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, EasyPlate, EasySep, RapidSpheres, and RoboSep are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.