

Small Molecules

Sodium Butyrate

Epigenetic modifier; Inhibits histone deacetylase

Catalog # 72242

500 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Sodium Butyrate is the sodium salt of butyric acid, a short chain fatty acid that inhibits histone deacetylases (HDACs), leading to hyperacetylation of histones. This causes changes in chromatin structure and gene expression, resulting in many biological effects (Boffa et al.; Kruh).

Molecular Name:	Sodium Butyrate
Alternative Names:	Butanoic acid; Butyric acid, sodium salt
CAS Number:	156-54-7
Chemical Formula:	$C_4H_7O_2 \cdot Na$
Molecular Weight:	110.1 g/mol
Purity:	≥ 95%
Chemical Name:	Not applicable
Structure:	

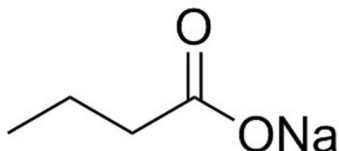

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at room temperature (15 - 25°C) as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.
Solubility:	· PBS (pH 7.2) ≤ 90 mM · Absolute ethanol ≤ 45 mM For example, to prepare a 10 mM stock solution in PBS, resuspend 100 mg in 90.8 mL of PBS (pH 7.2).

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in absolute ethanol at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final ethanol concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Supports self-renewal of mouse and human ES cells, in the absence of exogenously added growth factors (Ware et al.).

REPROGRAMMING

- Promotes reprogramming of human somatic cells to induced pluripotent stem (iPS) cells using only a single factor, OCT4 (Zhu et al.).

DIFFERENTIATION

- Promotes differentiation to hepatocytes from mouse and human embryonic stem (ES) cells (Hay et al., Zhou et al.).
- Promotes differentiation to definitive endoderm and islet-like cells from human ES cells (Jiang et al.).
- Enhances osteogenic and suppresses adipogenic differentiation from human mesenchymal cells (Chen et al., Lee et al.).

References

- Boffa LC et al. (1978) Suppression of histone deacetylation in vivo and in vitro by sodium butyrate. *J Biol Chem* 253(10): 3364–6.
- Chen T-H et al. (2007) Sodium butyrate activates ERK to regulate differentiation of mesenchymal stem cells. *Biochem Biophys Res Commun* 355(4): 913–8.
- Hay DC et al. (2008) Efficient differentiation of hepatocytes from human embryonic stem cells exhibiting markers recapitulating liver development in vivo. *Stem Cells* 26(4): 894–902.
- Jiang J et al. (2007) Generation of insulin-producing islet-like clusters from human embryonic stem cells. *Stem Cells* 25(8): 1940–53.
- Kruh J. (1982) Effects of sodium butyrate, a new pharmacological agent, on cells in culture. *Mol Cell Biochem* 42(2): 65–82.
- Lee S et al. (2009) Histone deacetylase inhibitors decrease proliferation potential and multilineage differentiation capability of human mesenchymal stem cells. *Cell Prolif* 42(6): 711–20.
- Ware CB et al. (2009) Histone deacetylase inhibition elicits an evolutionarily conserved self-renewal program in embryonic stem cells. *Cell Stem Cell* 4(4): 359–69.
- Zhou M et al. (2010) Differentiation of mouse embryonic stem cells into hepatocytes induced by a combination of cytokines and sodium butyrate. *J Cell Biochem* 109(3): 606–14.
- Zhu S et al. (2010) Reprogramming of human primary somatic cells by OCT4 and chemical compounds. *Cell Stem Cell* 7(6): 651–5.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2015 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.