

**Anti-Human CD184 (CXCR4)
Antibody, Clone 12G5, PE**



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Antibodies

Mouse monoclonal IgG2a antibody
against human, rhesus, cynomolgus
CD184 (CXCR4), PE-conjugated

Catalog #60089PE
#60089PE.1

100 Tests	5 µL/test
25 Tests	5 µL/test

Product Description

The 12G5 antibody reacts with CD184 (CXCR4 or fusin), an ~45 kDa seven transmembrane G-protein-coupled receptor expressed on the surface of a broad range of cells, including T and B cells, dendritic cells, monocytes, macrophages, granulocytes, platelets, endothelial and epithelial cells, astrocytes and neurons. At least two isoforms of CD184 are generated by alternate mRNA splicing. Ligands for CD184 include the CXC chemokine SDF-1 and ubiquitin. CD184 has functional roles in mediating chemotaxis, B cell lymphopoiesis and myelopoiesis, cardiogenesis, vascularization and neuronal guidance. CD184 acts as a co-receptor for infection by some T cell-tropic isolates of HIV-1 and as an alternative receptor for CD4-independent infection by HIV-2; functions that are blocked by binding of the 12G5 antibody. The antibody also partially inhibits SDF-1-induced chemotaxis and calcium influx. In differentiation studies, CD184 is used along with other markers such as CD117 to monitor the emergence of definitive endoderm.

Target Antigen Name:	CD184 (CXCR4)
Alternative Names:	C-X-C chemokine receptor type 4, CXCR4, fusin, HM89, HUMSTR, LAP3, LCR1, LESTR, NPY3R, NPY3R, SDF-1 receptor, WHIM
Gene ID:	7852
Species Reactivity:	Human, Rhesus, Cynomolgus, Baboon, Chimpanzee, African Green Monkey, Sooty Mangabey
Host Species:	Mouse (BALB/c)
Clonality:	Monoclonal
Clone:	12G5
Isotype:	IgG2a, kappa
Immunogen:	CP-MAC virus-infected SUP-T1 cells of human origin
Conjugate:	PE

Applications


Verified:	FC
Reported:	FC
Special Applications:	This antibody clone has been verified for purity assessments of cells cultured with STEMdiff™ Definitive Endoderm Kit (Catalog #05110).

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FACS: Fluorescence activated cell sorting; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; RIA: Radioimmunoassay; WB: Western blotting

Properties

Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) bovine serum albumin
Purification:	The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions. The solution is free of unconjugated PE and unconjugated antibody.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. Protect product from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com .
Directions for Use:	For flow cytometry the suggested use of this antibody is 5 µL per 1 x 10 ⁶ cells in 100 µL volume or per 100 µL of whole blood. It is recommended that the antibody be titrated for optimal performance for each application.

Data

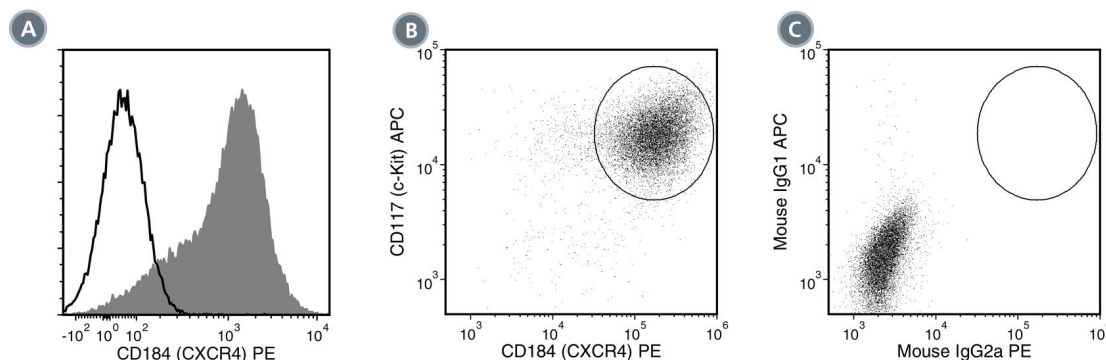

(A) Flow cytometry analysis of human peripheral blood mononuclear cells (PBMCs) labeled with Anti-Human CD184 (CXCR4) Antibody, Clone 12G5, PE (filled histogram) or Mouse IgG2a, kappa Isotype Control Antibody, Clone MOPC-173, PE (Catalog #60071PE) (solid line histogram).

(B) Flow cytometry analysis of human embryonic stem (ES) cells following differentiation to definitive endoderm using the STEMdiff™ Definitive Endoderm Kit. Cells were labeled with Anti-Human CD117 (c-Kit) Antibody, Clone 104D2, APC (Catalog #60087AZ) and Anti-Human CD184 (CXCR4) Antibody, Clone 12G5, PE. The emergence of definitive endoderm can be monitored by co-expression of CD117 and CD184.

(C) Flow cytometry analysis of human ES cells following differentiation to definitive endoderm using the STEMdiff™ Definitive Endoderm Kit. Cells were labeled with a mouse IgG1, kappa APC isotype control antibody and Mouse IgG2a, kappa Isotype Control Antibody, Clone MOPC-173, PE.

Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/antibodies or contact us at techsupport@stemcell.com.

References

1. Espirito-Santo M et al. (2012) Susceptibility of HIV type 2 primary isolates to CCR5 and CXCR4 monoclonal antibodies, ligands, and small molecule inhibitors. *AIDS Res Hum Retroviruses* 28(5): 478–85. (FA)
2. Uy GL et al. (2012) A phase 1/2 study of chemosensitization with the CXCR4 antagonist plerixafor in relapsed or refractory acute myeloid leukemia. *Blood* 119(17): 3917–24. (FC)
3. Grundle R et al. (2009) Dissection of PIM serine/threonine kinases in FLT3-ITD-induced leukemogenesis reveals PIM1 as regulator of CXCL12-CXCR4-mediated homing and migration. *J Exp Med* 206(9): 1957–70. (FC, ICC, IF)
4. Murga M et al. (2005) Neuropilin-1 regulates attachment in human endothelial cells independently of vascular endothelial growth factor receptor-2. *Blood* 105(5): 1992–9. (FC, ICC, IF)
5. Weng AP et al. (2003) CXCR4/CD184 immunoreactivity in T-cell non-Hodgkin lymphomas with an overall Th1- Th2+ immunophenotype. *Am J Clin Pathol* 119(3): 424–30. (IHC)
6. Yoshino N et al. (2000) Upgrading of flow cytometric analysis for absolute counts, cytokines and other antigenic molecules of cynomolgus monkeys (*Macaca fascicularis*) by using anti-human cross-reactive antibodies. *Exp Anim* 49(2): 97–110. (FC)
7. Berndt C et al. (1998) CXCR4 and CD4 mediate a rapid CD95-independent cell death in CD4(+) T cells. *Proc Natl Acad Sci USA* 95(21): 12556–61. (FA/Blocking)
8. Bleul CC et al. (1997) The HIV coreceptors CXCR4 and CCR5 are differentially expressed and regulated on human T lymphocytes. *Proc Natl Acad Sci USA* 94(5): 1925–30.
9. Endres MJ et al. (1996) CD4-independent infection by HIV-2 is mediated by fusin/CXCR4. *Cell* 87(4): 745–56. (IHC)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and STEMdiff are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.