

Antibodies

Mouse IgG3, kappa Isotype Control Antibody, Clone MG3-35

Mouse monoclonal IgG3, kappa
isotype control antibody,
unconjugated



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #60073
#60073.1

500 µg 0.5 mg/mL
50 µg 0.5 mg/mL

FOR RESEARCH USE ONLY. NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES.

Product Description

The MG3-35 antibody (IgG3, kappa) is suitable for use as an isotype-matched control antibody in several applications to estimate the degree of non-specific binding by an antigen-specific antibody. Ideally, the isotype control should have the same subclass of heavy chain (IgA, IgD, IgE, IgG, or IgM) and light chain (kappa or lambda) as the specific antibody being employed. If a conjugated antibody is employed, an isotype control conjugated to the same molecule (e.g. fluorochrome) should be chosen. The use of an appropriate isotype control helps confirm the specificity of the antigen-specific antibody and indicates non-specific binding that may result from binding to Fc receptors or other cell components. The MG3-35 antibody recognises keyhole limpet hemocyanin and has unknown binding specificity, having been screened on a variety of activated, resting, live, and fixed tissues from several species, including mouse, rat, human, and non-human primates.

Target Antigen Name:	IgG3 Isotype Control
Alternative Names:	Not applicable
Gene ID:	Not applicable
Species Reactivity:	Not applicable
Host Species:	Mouse
Clonality:	Monoclonal
Clone:	MG3-35
Isotype:	IgG3, kappa
Immunogen:	Trinitrophenol + KLH
Conjugate:	Unconjugated

Applications


Verified:	FC, ICC, IF
Reported:	FA, FC, ICC, IF, IHC, IP, WB
Special Applications:	This antibody clone has been verified for use as an isotype control antibody for assessing non-specific binding to cells in flow cytometry and immunofluorescence microscopy applications (surface and intracellular staining).

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; RIA: Radioimmunoassay; WB: Western blotting

Properties

Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Purification:	The antibody was purified by affinity chromatography.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. For product expiry date, please contact techsupport@stemcell.com .
Directions for Use:	The suggested use of this antibody is at concentrations comparable to those of the specific antibody of interest.

Data

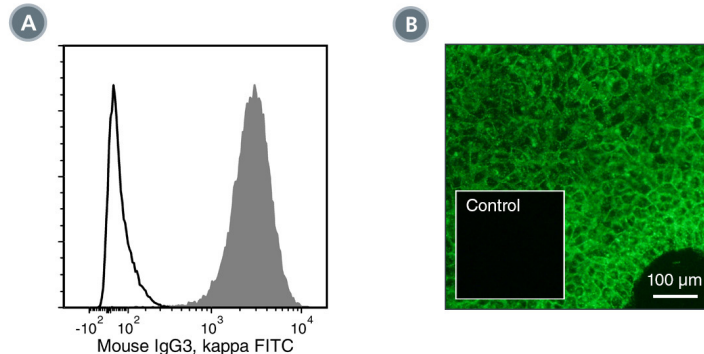

(A) Flow cytometry analysis of human iPS cells labeled with Mouse IgG3, kappa Isotype Control Antibody, Clone MG3-35, followed by a goat anti-mouse IgG antibody, FITC (solid line histogram). Filled histogram shows labeling with a mouse IgG3, kappa positive control antibody (Anti-Human SSEA-4 Antibody, Clone MC-813-70; Catalog #60062) followed by a goat anti-mouse IgG antibody, FITC.

(B) Human ES cells were cultured in mTeSR™1 on Corning® Matrigel®-coated glass slides, then fixed and stained with Anti-Human SSEA-4 Antibody, Clone MC-813-70 (Catalog #60062), followed by goat anti-mouse IgG, FITC. Inset shows cells labeled with Mouse IgG3, kappa Isotype Control Antibody, Clone MG3-35.

Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/antibodies or contact us at techsupport@stemcell.com.

References

1. Ameres S et al. (2014) CD8 T cell-evasive functions of human cytomegalovirus display pervasive MHC allele specificity, complementarity, and cooperativity. *J Immunol* 192(12): 5894–905. (FC)
2. Gao C et al. (2014) Carbohydrate sequence of the prostate cancer-associated antigen F77 assigned by a mucin O-glycome designer array. *J Biol Chem* 289(23): 16462–77. (FA)
3. De Riva A et al. (2013) Accelerated turnover of MHC class II molecules in nonobese diabetic mice is developmentally and environmentally regulated in vivo and dispensable for autoimmunity. *J Immunol* 190(12): 5961–71. (FC)
4. Fewou SN et al. (2012) Anti-ganglioside antibody internalization attenuates motor nerve terminal injury in a mouse model of acute motor axonal neuropathy. *J Clin Invest* 122(3): 1037–51. (FA)
5. Ribot JC et al. (2012) B7-CD28 costimulatory signals control the survival and proliferation of murine and human gamma delta T cells via IL-2 production. *J Immunol* 189(3): 1202–8. (FC)
6. Mitchell J & Sullam PM. (2009) Streptococcus mitis phage-encoded adhesins mediate attachment to α 2-8-linked sialic acid residues on platelet membrane gangliosides. *Infect Immun* 77(8): 3485–90. (FA)
7. Duan J et al. (2008) Microbial carbohydrate depolymerization by antigen-presenting cells: deamination prior to presentation by the MHC class II pathway. *Proc Natl Acad Sci USA* 105(13): 5183–8. (FA)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485 MEDICAL DEVICE STANDARDS.

Copyright © 2015 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and EasySep are trademarks of STEMCELL Technologies Inc. mTeSR is a trademark of WARF. Corning and Matrigel are trademarks of Corning Incorporated. All other trademarks are the property of their respective holders. Alexa Fluor® is a registered trademark of Life Technologies Corporation. This product is licensed for internal research use only and its sale is expressly conditioned on the buyer not using it for manufacturing, performing a service, or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.