

Cytokines

Human Recombinant Shh (C24II)

Sonic hedgehog



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #	78065.1	10 µg
	78065	50 µg
	78065.2	1000 µg

Product Description

Sonic hedgehog (Shh) is a member of the Hedgehog family of secreted signaling proteins that has a prominent role in patterning during early mammalian development (Ho & Scott). It is expressed during early embryogenesis in a variety of tissues (Finco et al.). During patterning of the central nervous system, Shh directly acts on cells to specify neural cell fates, regulates proliferation and survival of oligodendrocytes precursors and neural crest cells and stimulates proliferation, differentiation and axon growth (Ho & Scott; Marti). Shh also regulates T cell differentiation and activation (Crompton et al.). Recombinant human Shh (C24II; rhShh) is a fully biologically active molecule.

Product Information

Alternative Names:	HHG-1, HLP3, HPE3, MCOPCB5, SMMCI, Sonic hedgehog, TPT, TPTPS
Accession Number:	Q15465
Amino Acid Sequence:	IIGPGRGFGK RRHPKKLTPL AYKQFIPNVA EKT LGASGRY EGKISRNSER FKELTPNYPN DIIFKDEENT GADRLMTQRC KDKLNALAI VMNQWPGVKL RVTEGWDEDG HHSEESLHYE GRAVDITTSD RDRSKYGMLA RLAVEAGFDW VYYESKAHIH CSVKAENSVA AKSGG
Predicted Molecular Mass:	19.7 kDa
Species:	Human
Cross Reactivity:	Mouse
Formulation:	Lyophilized after dialysis against phosphate-buffered saline.
Source:	E. coli


Specifications

Activity:	The specific activity is ≥ 500 units/mg ($EC_{50} \leq 2$ µg/mL) as is determined by inducing alkaline phosphatase production of C3H/10T1/2 (CCL-226) cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic limulus amoebocyte lysate (LAL) analysis and is ≤ 0.2 EU/µg protein.

Preparation and Storage

Storage:	Store at -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Reconstitution:	Centrifuge vial before opening. Resuspend the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 2 weeks or at -20°C to -80°C for up to 3 months. Avoid repeated freeze-thaw cycles.

Data

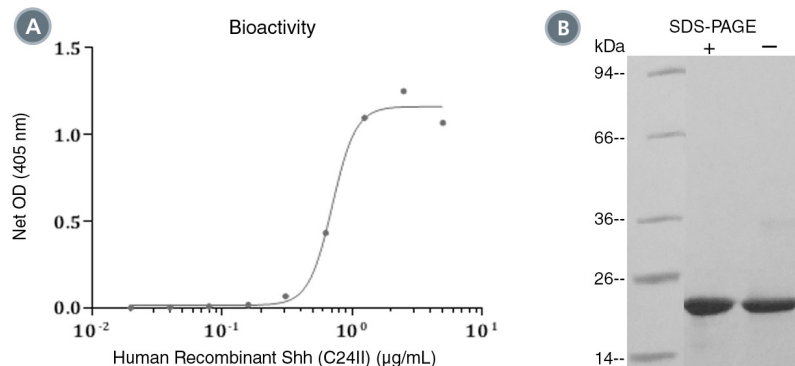

(A) The biological activity of Human Recombinant Shh (C24II) was tested by its ability to promote alkaline phosphatase production of CCL-226 cells. Alkaline phosphatase production was measured using a fluorometric assay method. The EC50 is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC50 in the above example 0.7 μg/mL.

(B) 4 μg of Human Recombinant Shh (C24II) was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant Shh (C24II) has a predicted molecular mass of 19.7 kDa.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Crompton T et al. (2007) Sonic hedgehog signalling in T-cell development and activation. *Nat Rev Immunol* 7(9): 726–35.
- Finco I et al. (2015) Hedgehog signaling and steroidogenesis. *Annu Rev Physiol* 77: 105–29.
- Ho KS & Scott MP. (2002) Sonic hedgehog in the nervous system: functions, modifications and mechanisms. *Curr Opin Neurobiol* 12(1): 57–63.
- Martí E. (2002) Sonic hedgehog in CNS development: one signal, multiple outputs. *Trends Neurosci* 25(2): 89–96.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.